

Untitled study

Glucose (mg/dL) · Specimen stability × Criteria not met

3 of 3 questions judged · Tested duration within the allowable change (time unit) and Claimed duration not met · Tested duration within the allowable change (time unit) met

Study date **Not recorded**

Table 1. Acceptance criteria

Question	Limit	Allowed	Finding	Result
Tested duration within the allowable change (time unit)	≥ 24 h		24 h	✓ Met
Tested duration within the allowable change (time unit)	≥ 48 h		24 h	× Not met
Claimed duration	≥ 48 h		24 h	× Not met ^a

^a Supported to 24 h. 48 h is beyond the allowable change.

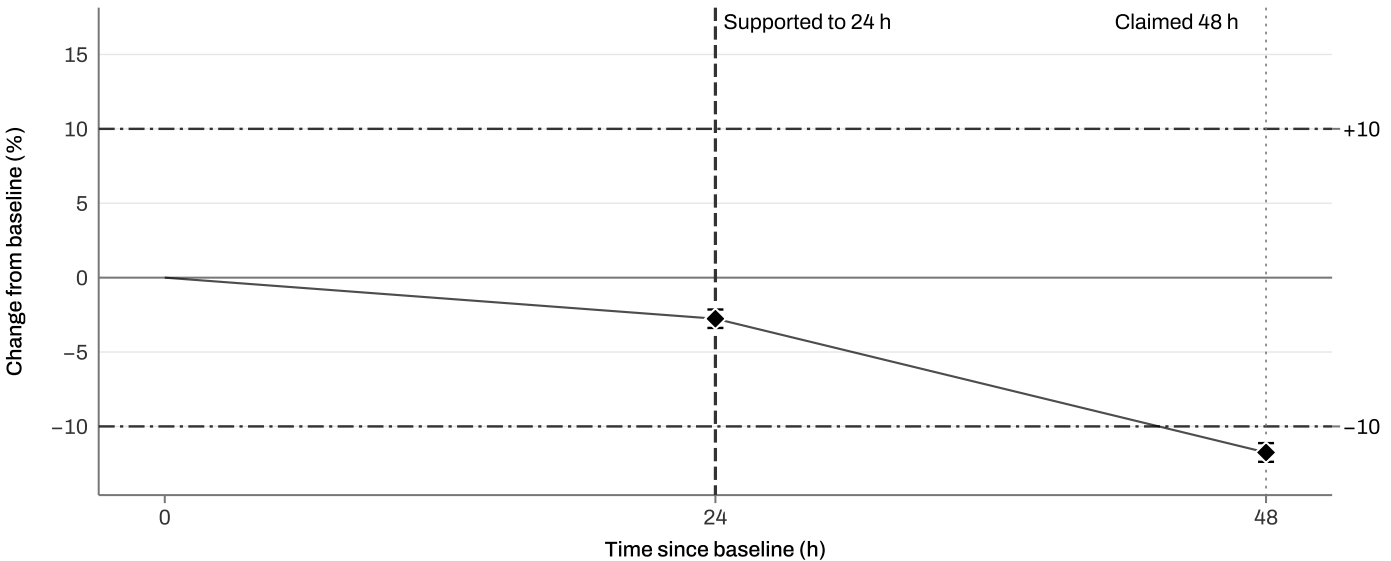


Figure 1. Ambient: mean change from baseline with 95% CI.

Table 2. Change from baseline by timepoint

Time (h)	N	Change (%)	95% CI (%)	Decision
24	5	-2.8	-3.4 to -2.1	✓ Met
48	5	-11.8	-12.4 to -11.1	× Not met

Study details

Data	30 included	Allowable change	±10%	Conditions	1 · 2 timepoints
------	-------------	------------------	------	------------	------------------

Signed off by

Date

A Data listing

Source: manual entry. 30 rows, all included.

Study row	Material	Condition	Time (h)	Result (mg/dL)	Replicate
1	S1	Shared baseline	0	82	1
2	S1	Shared baseline	0	82	2
3	S1	Ambient	24	81	1
4	S1	Ambient	24	79	2
5	S1	Ambient	48	73	1
6	S1	Ambient	48	73	2
7	S2	Shared baseline	0	95	1
8	S2	Shared baseline	0	95	2
9	S2	Ambient	24	93	1
10	S2	Ambient	24	93	2
11	S2	Ambient	48	84	1
12	S2	Ambient	48	84	2
13	S3	Shared baseline	0	119	1
14	S3	Shared baseline	0	118	2
15	S3	Ambient	24	114	1
16	S3	Ambient	24	115	2
17	S3	Ambient	48	104	1
18	S3	Ambient	48	104	2
19	S4	Shared baseline	0	140	1
20	S4	Shared baseline	0	140	2
21	S4	Ambient	24	135	1
22	S4	Ambient	24	137	2
23	S4	Ambient	48	123	1
24	S4	Ambient	48	123	2
25	S5	Shared baseline	0	164	1
26	S5	Shared baseline	0	166	2
27	S5	Ambient	24	160	1
28	S5	Ambient	24	160	2
29	S5	Ambient	48	145	1
30	S5	Ambient	48	146	2

1 Results

1.1 Timepoint detail

Table 3. Results behind each timepoint and the change in mg/dL

Time (h)	Results	Baseline results	Change (mg/dL)	95% CI (mg/dL)
Ambient				
24	10	10	-3.4	-5.1 to -1.7
48	10	10	-14.2	-19.5 to -8.9

1.2 Change by specimen

Table 4. Change by specimen, Ambient (%)

Specimen	Baseline (mg/dL)	24 h	48 h
S1	82.0	-2.4	-11.0
S2	95.0	-2.1	-11.6
S3	118.5	-3.4	-12.2
S4	140.0	-2.9	-12.1
S5	165.0	-3.0	-11.8

1.3 Conditions

Table 5. Baseline of each condition

Condition	Baseline	Specimens	Baseline mean (mg/dL)
Ambient	Shared	5	120.1